

Laboratory Protocol

# **Characterization of Protein Stability, Catalysis and Aggregation**

**Molecular Biology and Biochemistry of Life**  
Winter Semester 2023-24

**Authors**

Masine Janet Tränkner  
Syed Alisamar Husain

**Supervisor**

Dr. Titus M. Franzmann

## Contents

1	Temperature Induced Unfolding Transitions	2
2	Lysozyme Aggregation Assay	4
3	Michaelis-Menten Kinetics	5
4	UV/VIS Spectroscopy of NADH	7

## List of Figures

1	Lysozyme Fluorescence at wavelengths (a) 350nm and (b) 330nm . . . . .	2
2	BSA Fluorescence at wavelengths (a) 350nm and (b) 330nm . . . . .	2
3	Ratio of Fluorescence at 350nm to 330nm for Lysozyme . . . . .	3
4	Ratio of Fluorescence at 350nm to 330nm for BSA . . . . .	3
5	Scattering recorded as a function of time. . . . .	4
6	Rates . . . . .	5
7	Rates . . . . .	6
8	Absorbance spectrum of NADH for different concentrations . . . . .	7
9	Absorbance at 340nm for different concentrations . . . . .	7

## Author Information

Masine Janet Tränkner  
MSc. Computational Modelling and Simulation  
Matriculation Number 5186813  
*masine.janet.traenkner@mailbox.tu-dresden.de*

Syed Alisamar Husain  
MSc. Computational Modelling and Simulation  
Matriculation Number 5198172  
*syed.alisamar.husain@mailbox.tu-dresden.de*

# 1 Temperature Induced Unfolding Transitions

The goal in this experiment is to determine the stability of the proteins Lysozyme and BSA (Bovine Serum Albumin) at different concentrations, 2.5, 5 and 10  $\mu\text{M}$ , by inducing unfolding transitions through temperature changes. We detect unfolding using Fluorescence spectroscopy by measuring emission at two wavelengths, 330nm and 350nm, and we expect to see a sharp increase in the ratio of these emissions when the protein unfolds.

## Measurements

The temperature was increased from 25  $^{\circ}\text{C}$  to 95  $^{\circ}\text{C}$ , at a rate of 2  $^{\circ}\text{C}$  per minute, and corresponding fluorescence at wavelengths 350nm and 330nm was measured.

- **Lysozyme - Fluorescence vs. Temperature**

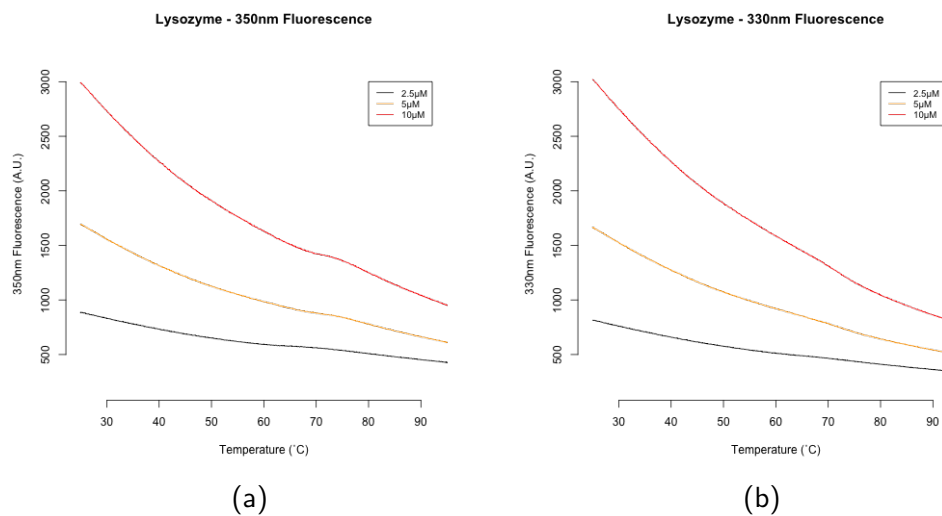


Figure 1: Lysozyme Fluorescence at wavelengths (a) 350nm and (b) 330nm

- **BSA - Fluorescence vs. Temperature**

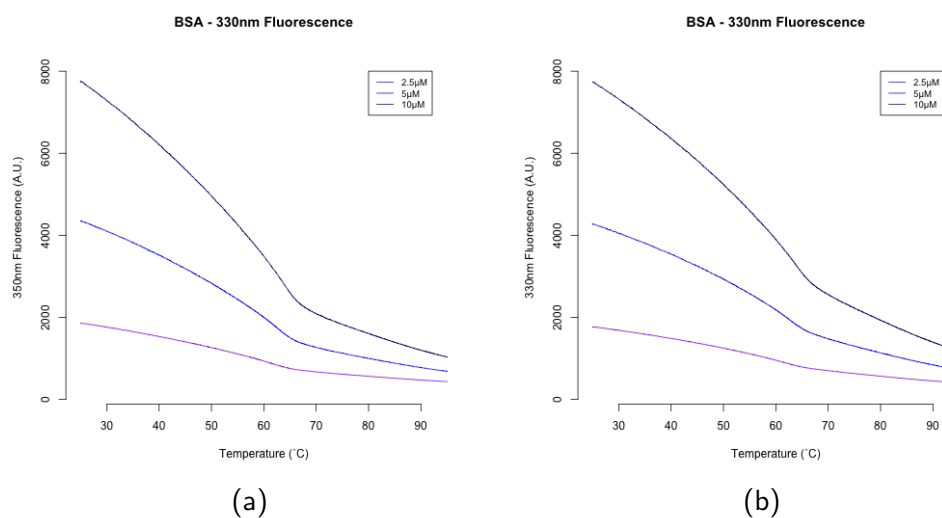


Figure 2: BSA Fluorescence at wavelengths (a) 350nm and (b) 330nm

## Analysis

The Fluorescence Ratio graph (Figure 3 and 4) was used to approximate the transition midpoint, while the protein's heat capacity was derived from the plot's derivative. Stability in the folded state arises from various interactions within the polypeptide. However, upon surpassing a threshold temperature, the polypeptide transitions to an unstable state. Consequently, as temperature increases, the population of unfolded molecules also rises. This process behaves differently for the two proteins.

For Lysozyme, a low ratio can be observed at the beginning, which indicates that the folded state predominates. With increasing temperature, this confirmation changes up to a threshold temperature and the transition of the molecules into the unfolded state, and accordingly the ratio also increases.

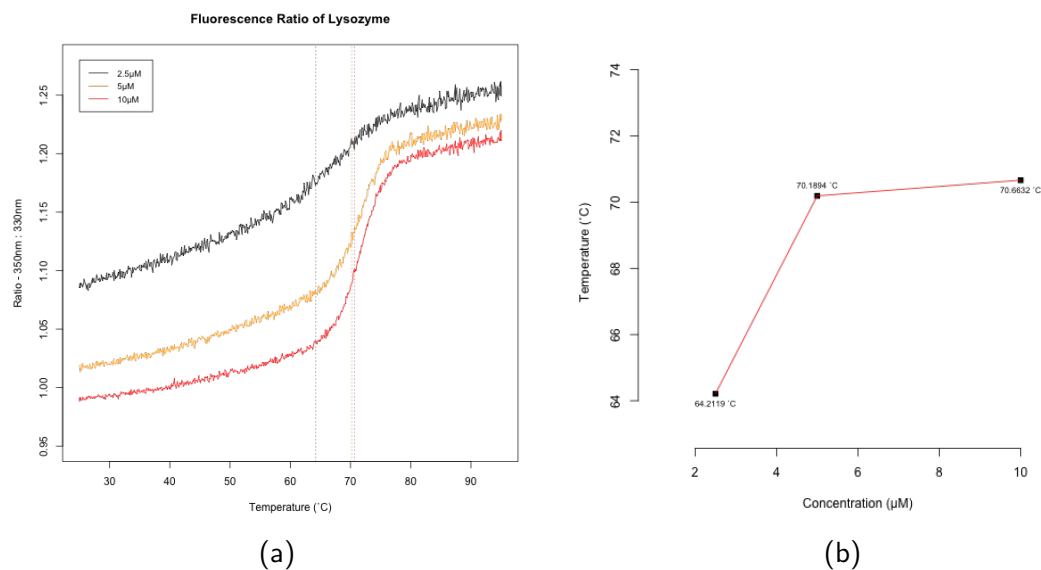


Figure 3: Ratio of Fluorescence at 350nm to 330nm for Lysozyme

The opposite can be observed for BSA. Here, the ratio is already high at the start of the process and falls down to a certain threshold temperature. From this point, the ratio increases for all concentrations and thus also the amount of proteins in the unfolded state.

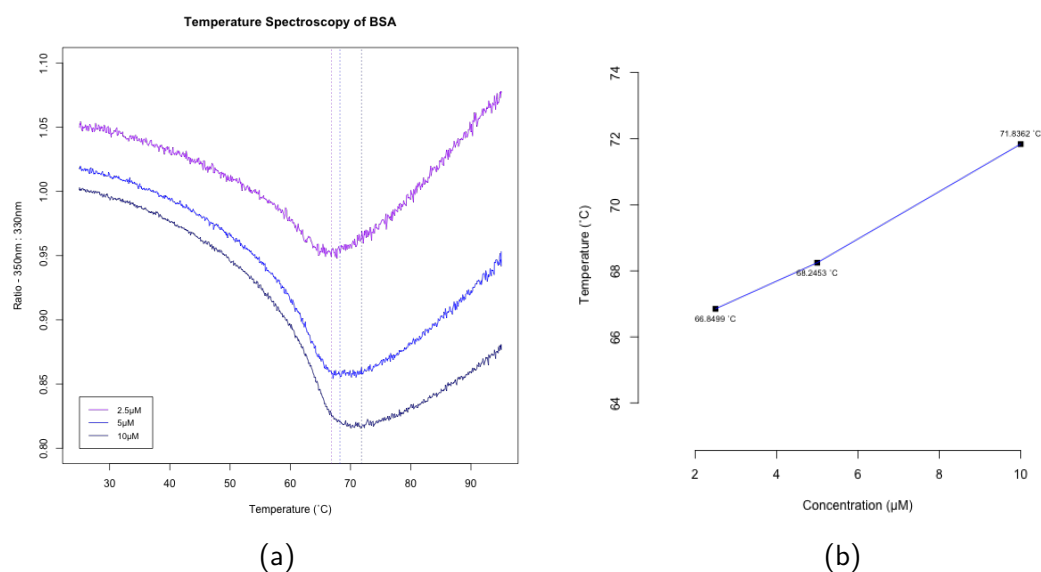


Figure 4: Ratio of Fluorescence at 350nm to 330nm for BSA

The transition midpoints are plotted in the two plots.

## 2 Lysozyme Aggregation Assay

In this experiment, aggregation of Lysozyme was studied to understand protein aggregation, where unfolded polypeptide chains of the proteins Lysozyme and Hsp26 associate irreversibly, leading to the formation of large particles that scatter light. Throughout the experiment, Hsp26 was used to analyze its effect on the aggregation of Lysozyme with presence of TCEP.

### Measurements

Poor absorbance spectrometer used to do a poor-man's version of scattering. The scattering is not actually detected. What we do is choose a wavelength where we know nothing in the mix absorbs anything and we use that fact to say all of the intensity reduction is because of the scattered light.

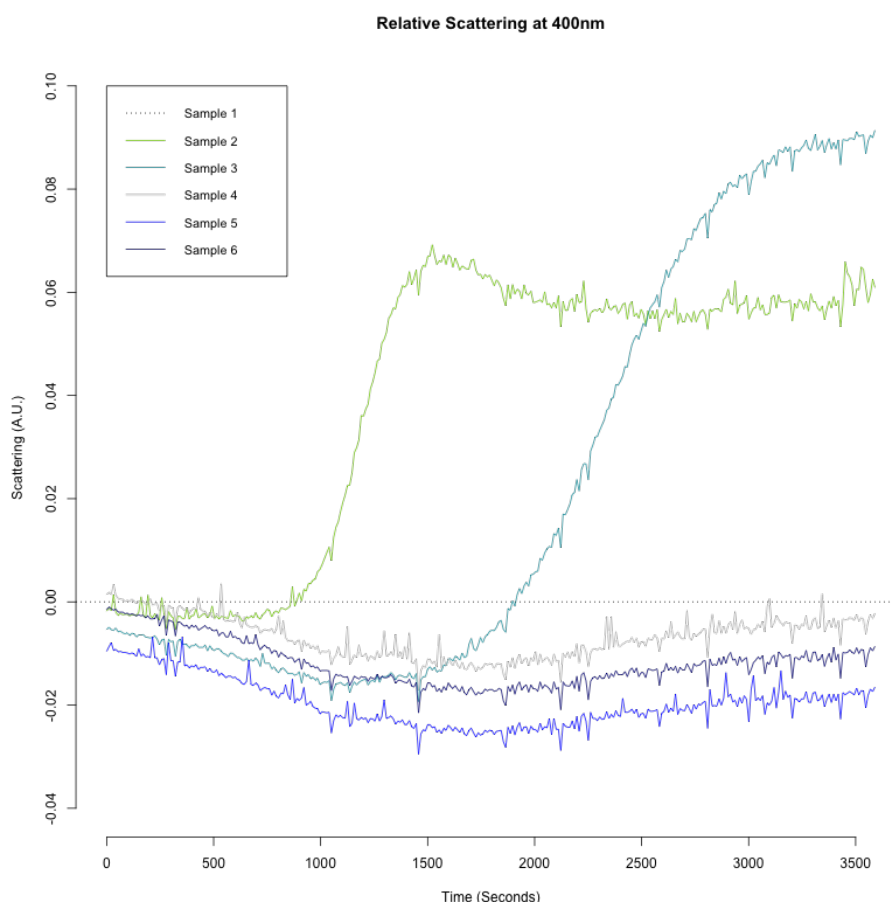


Figure 5: Scattering recorded as a function of time.

### Analysis

As protein aggregation increases, more light is absorbed by the sample. Introducing Hsp26 into the sample solutions can regulate aggregation and, consequently, absorption. At low Hsp26 concentrations (Samples 1 to 3 in Figure 5), Lysozyme aggregation isn't effectively inhibited. Thus, the presence of Hsp26 in the solution increases overall absorption. But it can also be seen that a higher concentration of 5  $\mu\text{M}$  delays the aggregation.

However, when the concentration of Hsp26 matches that of Lysozyme, aggregation is prevented, resulting in lower absorbance (Sample 4). The control samples lacking Lysozyme and denaturant exhibit no absorbance change at 400nm (Samples 5 and 6). Only Samples 2 and 3 allow for the determination of an apparent halftime, as their curves display a change of confirmation.

### 3 Michaelis-Menten Kinetics

Lactate dehydrogenase (LDH) catalyzes the reaction of the substrate pyruvate into the product, lactate. In order to measure its activity, the conversion of NADH, the cofactor of LDH,  $\text{NAD}^+$  was observed in this experiment through spectroscopy of NADH.

#### Measurements

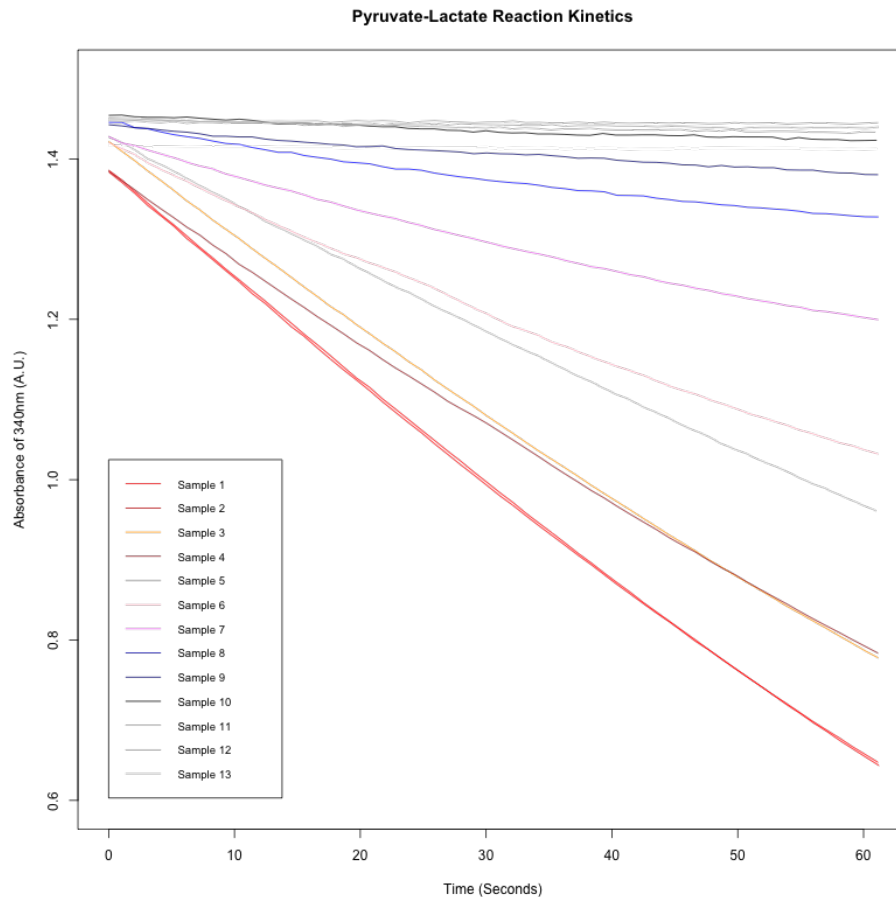


Figure 6: Rates

#### Analysis

The Michaelis-Menten plot (Figure 7(a)) could be plotted from the reaction rates and the according concentration of substrate (NADH). It shows that for low substrate concentrations, the rate of the reaction increases linearly with increasing substrate concentration. As the substrate concentration increases, the rate of the reaction or velocity eventually reaches a maximum value ( $V_{max}$ ). In this region, the active sites of the enzyme are largely occupied by the substrate, and the rate becomes independent of further increases in substrate concentration.

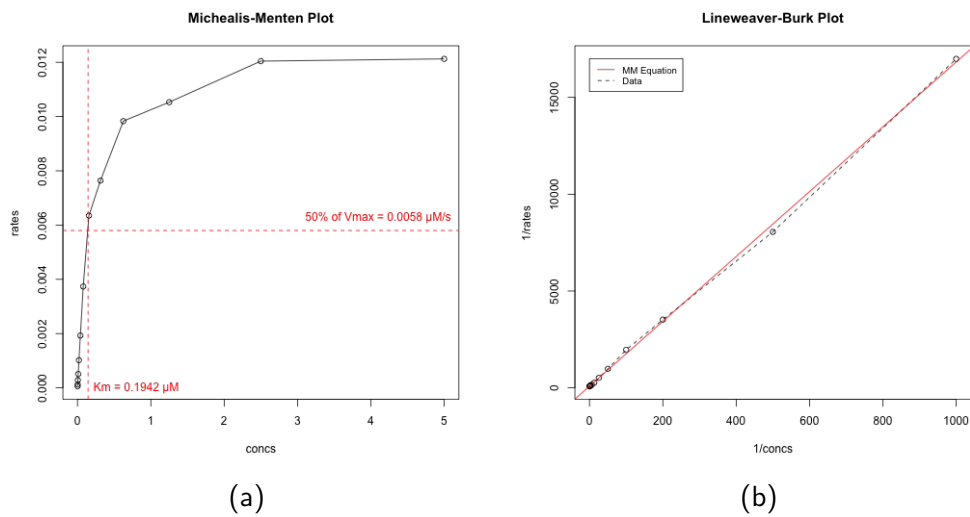


Figure 7: Rates

Derived from the Michaelis-Menten plot, the Lineweaver-Burk plot provides a linear transformation. It is used for determining the kinetic constants, the maximum reaction rate  $V_{max}$  and the Michaelis constant  $K_m$ .

The result for  $V_{max} = 0.0116 \mu\text{M/s}$  and for  $K_m = 0.1942 \mu\text{M}$ , respectively.

[S] $\mu\text{M}$	$\Delta\text{NADH}$ ( $\mu\text{M}$ )
5000	asf
2500	asf
1250	asf
5000	asf
5000	asf

## 4 UV/VIS Spectroscopy of NADH

Similar to the kinetics experiment, the concentration should also be measured here using UV/VIS spectroscopy. The detected absorbance should show local maxima for 3 specific wavelengths (208, 260 and 340 nm). It is also to be expected that the absorbance also increases with increasing NADH concentration. However, according to figure x, the absorbance increases from sample 1 to 3, then increases for the samples 4 and 5 and again increases for sample 6, without exceeding the highest absorbance of sample 2.

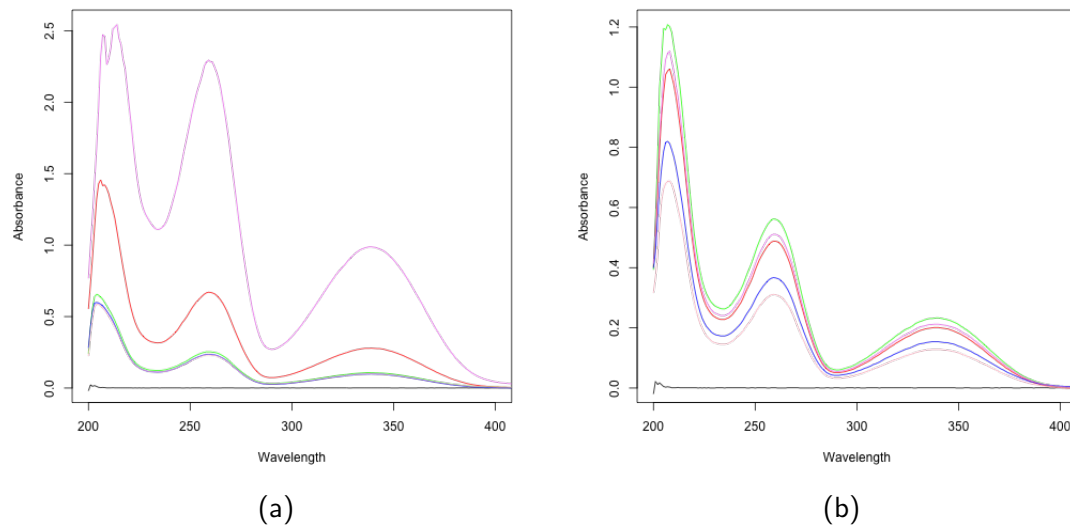


Figure 8: Absorbance spectrum of NADH for different concentrations

## Analysis

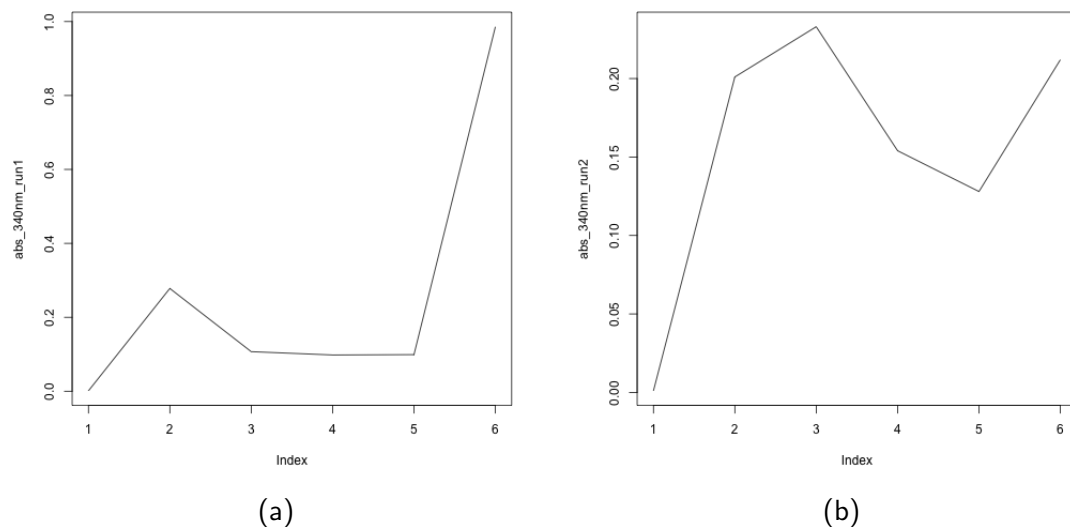


Figure 9: Absorbance at 340nm for different concentrations